

# Efficacy of Different C<sub>18</sub> HPLC Analytical Columns in the Analysis of Fumonisin B<sub>1</sub> and B<sub>2</sub> in Different Matrices

Yousef Sultan<sup>1,\*</sup> , Naresh Magan<sup>2</sup> , Angel Medina<sup>2</sup> 

<sup>1</sup> Food Toxicology and Contaminants Department, National Research Centre, Dokki, Cairo 12622, Egypt

<sup>2</sup> Applied Mycology Group, Cranfield Soil and AgriFood Institute, Cranfield University, Cranfield, Bedford MK43 0AL, United Kingdom

\* Correspondence: [sultanay18@hotmail.com](mailto:sultanay18@hotmail.com) (Y.S.);

Scopus Author ID 36480540200

Received: 11.04.2021; Revised: 15.05.2021; Accepted: 18.05.2021; Published: 10.06.2021

**Abstract:** Fumonisin B<sub>1</sub> and B<sub>2</sub> are carcinogenic and commonly contaminate corn and corn-based products. Analysis of such toxins using C<sub>18</sub> HPLC column is officially accredited but still unknown if all column types can effectively separate FB<sub>1</sub> and FB<sub>2</sub> or not. The present study evaluated the efficiency of 5 analytical columns with different dimensions, particle sizes, and porosities to determine these toxins in both agar cultures of *Fusarium verticillioides* and cornflakes. Interestingly, the traditional column 150mm of length with 5µm porous particles had close retention times to those of the short-fused core column 75mm of length with 2.7 µm reflecting in time and solvents saving. Using Sep-Pack C<sub>18</sub> for clean-up played an important role in enhancement the limit of quantification (LOQ) for cornflake samples (5-13.7 and 16.1-39 µg kg<sup>-1</sup> for FB<sub>1</sub> and FB<sub>2</sub>, respectively). However, it was relatively higher for fungal culture samples that were not passed through the cleaning-up step (11.5-16 and 28.1-46.3 µg kg<sup>-1</sup> for FB<sub>1</sub> and FB<sub>2</sub>, respectively). Overall, the lowest LOQ was obtained using the shorter fused core column. Finally, using such clean-up in the extraction of FB<sub>1</sub> and FB<sub>2</sub> from spiked cornflakes samples gave good recoveries (>80%) using all tested columns.

**Keywords:** fumonisin; particle size; column length; solid core; cornflakes.


© 2021 by the authors. This article is an open-access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (<https://creativecommons.org/licenses/by/4.0/>).

## 1. Introduction

Fumonisin are mycotoxins produced mainly by *Fusarium verticillioides* (formerly: *F. moniliforme*) and *Fusarium proliferatum* and are considered a global serious problem in corn crops [1-5].

Fumonisin B<sub>1</sub> (FB<sub>1</sub>) is 2-amino-12, 16-dimethyl-3,5,10,14,15-pentahydroxyeicosane esterified at C-14 and C-15 to propane-tricarboxylic acid, and fumonisin B<sub>2</sub> (FB<sub>2</sub>) is 10-deoxy-FB<sub>1</sub> (Figure 1, [1]). FB<sub>1</sub> has been classified as a group 2B carcinogen (possibly carcinogenic in humans) [6, 7]. It is found to cause equine leukoencephalomalacia and porcine pulmonary edema and a non-genotoxic kidney or liver cancer in rats and mice [8-11]. Several studies noticed a relation between FB<sub>1</sub> and human esophageal cancer [12, 13].

Fumonisin are mainly found as natural contaminants of corn and corn-based foods [1, 14-16]. However, there is evidence that they can occur in other crops and derived foods: sorghum [17] and sorghum syrup [18]; white beans, adzuki beans, and mung beans [19]; wheat, barley, and soybean [20, 21]; black tea and medicinal plants [22, 23]; rice [24].

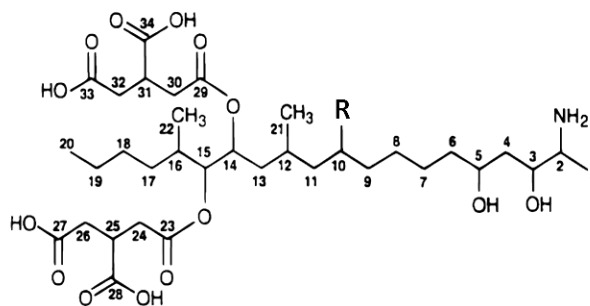

**Figure 1.** Structures of FB<sub>1</sub> (R=OH) and FB<sub>2</sub> (R=H).

Finding the appropriate method for fumonisins determination in foodstuff is the first step to establish a good control strategy since the severity of the occurrence needs to be determined. Such a method should also meet the international acceptance criteria limits [25].

The European Commission set up permissible limits of 4000  $\mu\text{g kg}^{-1}$  for unprocessed maize, 1000  $\mu\text{g kg}^{-1}$  for maize intended for direct human consumption, 800  $\mu\text{g kg}^{-1}$  in snacks & breakfast cereals, and 200  $\mu\text{g kg}^{-1}$  processed maize-based foods and babies and young children foods [26].

Based on its hydrophobicity properties, reversed-phase HPLC columns (RP, C<sub>18</sub>) are used to analyze most of the common toxins belonging to the fumonisins group [27]. Although several methods have been developed, either using mass spectrometry [1,28-32] or diode array [33] detectors, the fluorescence detector (FLD) is still the most common and appropriate for most matrices and still is, the official recommended method [34-36]. In addition, the analysis cost using HPLC equipped with FLD is cheaper than that using LC/MS, which is not available in many laboratories around the world.

Before the analysis, different clean-up techniques have been used to extract and purify fumonisins B<sub>1</sub>, B<sub>2</sub> from food samples. Most of these approaches have used strong anion exchange, SAX [14, 33], immunoaffinity column, IAC [37], or RP-C<sub>18</sub> solid-phase extraction (SPE) [20, 34, 38]. The lack of a native fluorophore in the fumonisin structure makes derivatization prior to the analysis necessary. Currently, the most common derivatizing agent used is o-phthalaldehyde (OPA) [36].

On the other hand, using new C<sub>18</sub> packing materials with fused solid core particles was found to improve the separation of aflatoxins [39] and ochratoxin A [40] by reducing the time and solvents used in the analysis with no effect on the performance.

Thus, this study aimed to examine the efficiency of 5 analytical columns with different particle sizes and lengths and C<sub>18</sub> chemistry to analyze FB<sub>1</sub> and FB<sub>2</sub> in synthetic culture media and natural and spiked cornflakes samples. The limit of detection (LOD), the limit of quantification (LOQ), and other chromatographic parameters (Peak width (W), number of theoretical plates (N), height equivalent to a theoretical plate (HEPT), and reduced plate height (h)) were used to perform the evaluation.

## 2. Materials and Methods

### 2.1. Reagents and standards.

A Fumonisins mixture (FB<sub>1</sub> and FB<sub>2</sub>, 50  $\mu\text{g ml}^{-1}$  Acetonitrile (CAN): water (1:1, v/v), Orthophthalaldehyde (OPA), orthophosphoric acid, sodium tetraborate (Na<sub>2</sub>B<sub>4</sub>O<sub>7</sub>), and 2-mercaptoethanol were purchased from Sigma-Aldrich, UK. Acetonitrile (AcN) and methanol (MeOH) were purchased from Fisher Scientific (Fisher Scientific UK Ltd., UK). C<sub>18</sub> Sep-Pak

solid-phase extraction (SPE) cartridges, 500 mg, were purchased from Variant (Union, Missouri, USA). All solvents were HPLC grade. The water used was obtained from a Milli-R/Q water system (Millipore, Billerica, MA, USA).

## 2.2. Working standards.

Calibration curves of FB<sub>1</sub> and FB<sub>2</sub> were prepared using a range of concentrations from 0.05 to 10 µg ml<sup>-1</sup> ACN:H<sub>2</sub>O, 1:1.

## 2.3. Media and sample preparation.

The medium used in this study was a fumonisin-inducing solid agar medium (FIM) previously used [41]. It contained for each liter 0.5 g malt extract, 1 g yeast extract, 1 g peptone, 1 g KH<sub>2</sub>PO<sub>4</sub>, 0.3 g MgSO<sub>4</sub>·7H<sub>2</sub>O, 0.3 g KCl, 0.05 g ZnSO<sub>4</sub>·7H<sub>2</sub>O, 0.01 g CuSO<sub>4</sub>·5H<sub>2</sub>O, 20 g fructose, and 15 g bacteriological agar. The culture medium was autoclaved for 20 minutes at 121°C, vigorously shaken, and poured into 9 cm diameter sterile Petri dishes.

FIM plates were centrally inoculated with 3 µl spore suspensions (1x10<sup>6</sup> spores ml<sup>-1</sup>) of two strains of *F. verticillioides*, S30 and S40 (kindly provided by Dr. Sejakoshi Mohale, Cranfield). Four replicates of each strain were incubated in the dark at 25°C for 10 days. Afterward, 5 agar discs, including the fungus and agar, were removed from these cultures using a cork borer (0.8 cm), transferred to a pre-weighed 2 ml Eppendorf tube, weighed again, and frozen at -20°C until further fumonisins analysis by HPLC-FLD.

Cornflakes samples (500g package, Kellogg's, UK) were purchased from a retail shop in the UK, mixed well, and split into two subsamples. The first subsample represented the control samples, and the second was spiked using the fumonisin B<sub>1</sub> and B<sub>2</sub> standards at different concentrations (0.625, 1.25, and 2.5 mg kg<sup>-1</sup>). Each subsample was subsequently divided into smaller samples (25g) for later analysis.

## 2.4. FB<sub>1</sub> and FB<sub>2</sub> extraction.

Toxins extraction from media was performed according to the method of Lazzaro *et al.* [42] by adding 1ml acetonitrile/water (1:1 v/v) to the plugs, shaking for 1h, and centrifugation for 5 min. The supernatant was subsequently filtered through a Nylon Filter 13 mm 0.22µm (Jaytee) in a new Eppendorf in preparation for derivatization.

Cornflakes samples were extracted according to Dombrink-Kurtzman and Dvorak [34]. Twenty-five grams of ground cornflakes were extracted with 100 ml of acetonitrile/water (1:1 v/v), shaken for 1 h, and filtered. A C<sub>18</sub> cartridge was preconditioned by passing through 5 ml of methanol and then 5 ml of water at 1-3 ml min<sup>-1</sup> flow rate. Two milliliters of sample filtrate were diluted with 5 ml of water and applied to the cartridge Sep-Pak. It was washed with 5 mL of water, followed by 2 ml of acetonitrile/water (1:9 v/v). The fumonisins were eluted with 4 ml of acetonitrile/water (7:3 v/v) at a flow rate of <1 ml min<sup>-1</sup>. The eluent was evaporated to dryness (under nitrogen), dissolved in 1 ml acetonitrile/water (1:1 v/v), filtered, and derivatized with OPA as described by AOAC [36].

## 2.5. Derivatization and LC analysis.

Fifty microliters aliquots of the sample extracts or standards were transferred to 250µl vial inserts (Agilent, Berks., UK), and 100µl of OPA reagent (40 mg OPA+1 ml methanol for dissolving+5ml 0.1M Na<sub>2</sub>B<sub>4</sub>O<sub>7</sub> +50µl 2-mercaptoethanol) were added and mixed by pipetting and injected to the HPLC system within 1 min.

## 2.6. Chromatographic equipment and fumonisin analysis.

The HPLC system used for FB<sub>1</sub>, FB<sub>2</sub> analysis was an Agilent 1200 Series system (Agilent, Berks., UK) equipped with a fluorescence detector (FLD, G1321A, Agilent), an autosampler (ALS, G1329, Agilent), autosampler thermostat (G1330B, Agilent), thermostatted column compartment (G1316A, Agilent), on-line degasser (G1379B, Agilent), and binary pump (G1312A, Agilent). The separations were performed in the isocratic mode. A slight modification of the mobile phase recommended by the AOAC official method [36] was used and contained 25% 0.1M NaH<sub>2</sub>PO<sub>4</sub>:75% methanol adjusted to pH 3.35 using phosphoric acid. The flow rate was set at 1 ml min<sup>-1</sup>. FLD detection was performed at 335nm excitation and 440nm emission wavelengths.

With an appropriate pre-column, 5 different columns were examined for their performance in FB<sub>1</sub> and FB<sub>2</sub> analysis. Three of them contained fully porous particles [Cronus Nucleosil 100 C<sub>18</sub> 4.6x150mm, 5µm; Agilent Zorbax Eclipse plus C<sub>18</sub> 4.6x150mm, 3.5µm and Agilent Zorbax Eclipse plus C<sub>18</sub> 4.6x100mm, 3.5µm], whereas Agilent Poroshell 120 EC-C<sub>18</sub> 4.6x100mm, 2.7µm and Agilent Poroshell 120 SB-C<sub>18</sub> 4.6x75mm, 2.7µm were packed with fused-core silica particles.

The injection volume was set at 20 µl. Agilent ChemStation software Ver. B Rev: 03.01 (Agilent Technologies, Palo Alto, CA, USA) was used to analyze the output signals.

## 2.7. Column efficiency.

The efficiency of all tested columns was assessed. Along with different runs, the backpressure of each column was recorded. In triplicate injections, the linearity of standard solutions was measured at concentrations of 0.66, 6.6, 66, and 134 ng of either FB<sub>1</sub> or FB<sub>2</sub> injection<sup>-1</sup>. Calibration lines were generated by plotting the peak area against the concentrations. Linear regression was obtained using Microsoft®Excel® in order to establish the correlation coefficient.

Based on the peak width (W) of either FB<sub>1</sub> or FB<sub>2</sub> obtained using Agilent ChemStation software, parameters of each column were calculated as such equations:

1)  $N = 16 \left( \frac{V_e}{W_b} \right)^2$ , where N is the number of theoretical plates, V<sub>e</sub> is the elution volume, and W<sub>b</sub> is the width of the peak at the baseline.

2)  $HEPT = L/N$ , where HEPT is the height equivalent theoretical plate and L is the column length

3)  $h = HEPT/dp$ , where h is the reduced plate height (h), and dp is the particle size (µm). That parameter is dimensionless, which facilitates the comparison of different columns packed with different particle sizes [43].

The sensitivity of the tested columns was assessed by the calculation of LOD and LOQ of FB<sub>1</sub> and FB<sub>2</sub> for standards, fungal cultures, and cornflake samples. LOD and LOQ were considered the fumonisins concentration that provided S/N ratio 3:1 and S/N ratio 10:1, respectively [44].

## 2.8. Statistical analyses.

Statistica Version 10 (StateSoft, Tulsa, Okla., USA) was used for data analysis. Analysis of variance (ANOVA, one way at  $p < 0.05$ ) was applied to fumonisin B<sub>1</sub> and B<sub>2</sub> concentrations in fungal cultures and cornflakes samples. Fisher's LSD method was used for the comparison between columns' parameters.

### 3. Results and Discussion

#### 3.1. The enhancement of FB1 and FB2 separation.

The first improvements were made using the Cronus 150x4.6mm-5 $\mu$ m column, and then the modified mobile phase was used with the other columns for the comparison in FB<sub>1</sub> and FB<sub>2</sub> separation (Figure 2).


FB<sub>1</sub> was totally separated, and no interference with any close peaks was noticed. The peak that appeared closer to FB<sub>2</sub> was also completely separated, lifting a completely resolute FB<sub>2</sub> peak.

#### 3.2. Inter-column comparison of the different chromatographic parameters.

Table 1 shows the data of column backpressure, total time of analysis, the retention time of both FB<sub>1</sub> and FB<sub>2</sub>, and the most important chromatographic parameters of the tested columns.

In general, the backpressures of all tested columns were within the safe range advised for most Agilent HPLC 1200 series (<400 Bar). An increase in the backpressure was noticed by decreasing the particle size of the same length columns (150mm: 5 and 3.5  $\mu$ m). Contrarily, a decrease in backpressure was observed by decreasing the length of the same particle size columns (100 and 75 mm: 2.7 $\mu$ m).

Similarly, the shorter the column was, the less the retention times obtained at the same particle sizes (150, 100mm: 3.5 $\mu$ m -100, 75mm: 2.7 $\mu$ m). At the length of 150 mm columns, the reduction in porous particle size from 5 to 3.5  $\mu$ m led to increased retention time (+3 and 11 min for FB<sub>1</sub> and FB<sub>2</sub>, respectively). However, the reduction from 3.5 $\mu$ m porous particles to 2.7  $\mu$ m solid core particles in columns of 100mm slightly decreased the retention times from 5.5 and 14.4 to 5 and 13.2 min for FB<sub>1</sub> and FB<sub>2</sub>, respectively.

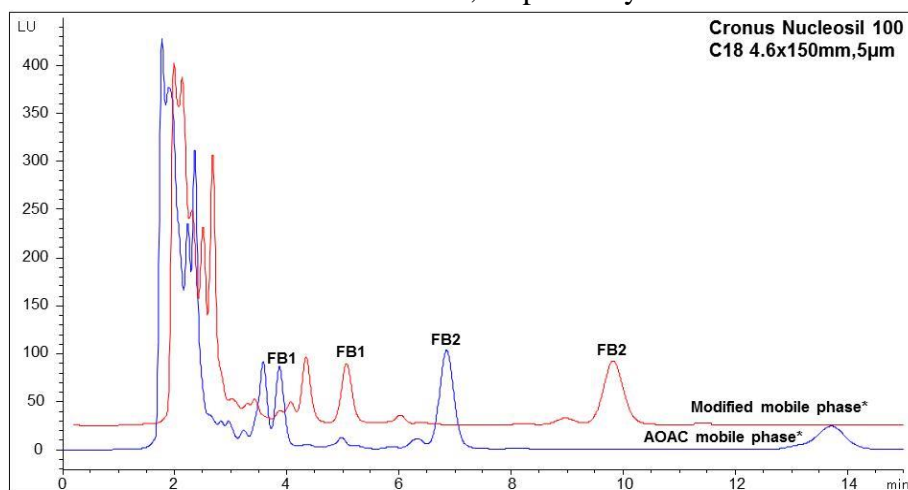

**Figure 2.** Enhancement of Fumonisin B<sub>1</sub> and B<sub>2</sub> separation in *Fusarium* culture medium using modified mobile phase.


**Table 1.** Backpressures, Fumonisin B<sub>1</sub> and B<sub>2</sub> (FB<sub>1</sub>, FB<sub>2</sub>) retention times ( $t_r$ ), and chromatographic parameters calculated for the different columns. Values obtained from 3 different injections at low, medium, and high concentrations. SE: Standard Error.

| C <sub>18</sub> Column     | Column pressure (Bar) | Toxin type      | $t_r \pm SE$      | $W \pm SE$          | $N \pm SE$           | HEPT $\pm SE$    | $h \pm SE$      | Total analysis time (min.) |
|----------------------------|-----------------------|-----------------|-------------------|---------------------|----------------------|------------------|-----------------|----------------------------|
| Cronus 150x4.6mm-5 $\mu$ m | 180 $\pm 2$           | FB <sub>1</sub> | 4.759 $\pm 0.002$ | 0.2068 $\pm 0.0038$ | 8488.16 $\pm 316.23$ | 17.72 $\pm 0.67$ | 3.54 $\pm 0.03$ | Cornflakes 105             |
|                            |                       | FB <sub>2</sub> | 9.269 $\pm 0.007$ | 0.3826 $\pm 0.0041$ | 9395.49 $\pm 214.47$ | 15.98 $\pm 0.37$ | 3.20 $\pm 0.07$ | Culture 90                 |

| C <sub>18</sub> Column         | Column pressure (Bar) | Toxin type      | t <sub>r</sub> ±SE | W±SE          | N±SE            | HEPT±SE   | h±SE       | Total analysis time (min.) |     |
|--------------------------------|-----------------------|-----------------|--------------------|---------------|-----------------|-----------|------------|----------------------------|-----|
| Zorbax Eclipse 150x4.6mm-3.5µm | 190±2                 | FB <sub>1</sub> | 7.791±0.006        | 0.2313±0.0034 | 18173.18±557.00 | 8.27±0.25 | 2.36±0.07  | Cornflakes                 | 115 |
|                                |                       | FB <sub>2</sub> | 20.649±0.019       | 0.5328±0.0087 | 24066.30±760.83 | 6.24±0.19 | 1.78±0.055 | Culture                    | 100 |
| Zorbax Eclipse 100x4.6mm-3.5µm | 132±2                 | FB <sub>1</sub> | 5.469±0.018        | 0.1988±0.0024 | 12126.33±374.39 | 8.26±0.25 | 2.36±0.07  | Cornflakes                 | 110 |
|                                |                       | FB <sub>2</sub> | 14.350±0.020       | 0.4329±0.0026 | 17587.01±249.09 | 5.69±0.08 | 1.63±0.02  | Culture                    | 95  |
| Poroshell 100x4.6mm-2.7µm      | 272±2                 | FB <sub>1</sub> | 4.993±0.006        | 0.1581±0.0019 | 15937.37±406.72 | 6.28±0.16 | 2.33±0.06  | Cornflakes                 | 110 |
|                                |                       | FB <sub>2</sub> | 13.156±0.016       | 0.3476±0.0007 | 22923.11±45.76  | 4.36±0.01 | 1.62±0.00  | Culture                    | 95  |
| Poroshell 75x4.6mm-2.7µm       | 232±2                 | FB <sub>1</sub> | 3.234±0.004        | 0.1241±0.0004 | 10825.42±58.83  | 6.93±0.04 | 2.57±0.01  | Cornflakes                 | 105 |
|                                |                       | FB <sub>2</sub> | 7.662±0.017        | 0.2402±0.0028 | 16295.46±397.59 | 4.61±0.11 | 1.71±0.04  | Culture                    | 90  |

W: peak width (min); N: number of theoretical plates; HEPT: height equivalent to a theoretical plate (µm); h: reduced plate height.

Figure 3 shows examples of the FLD chromatograms obtained by injecting both FB<sub>1</sub> and FB<sub>2</sub> standards, fungal culture extracts, and spiked and unspiked (control) cornflakes using the set of columns tested.

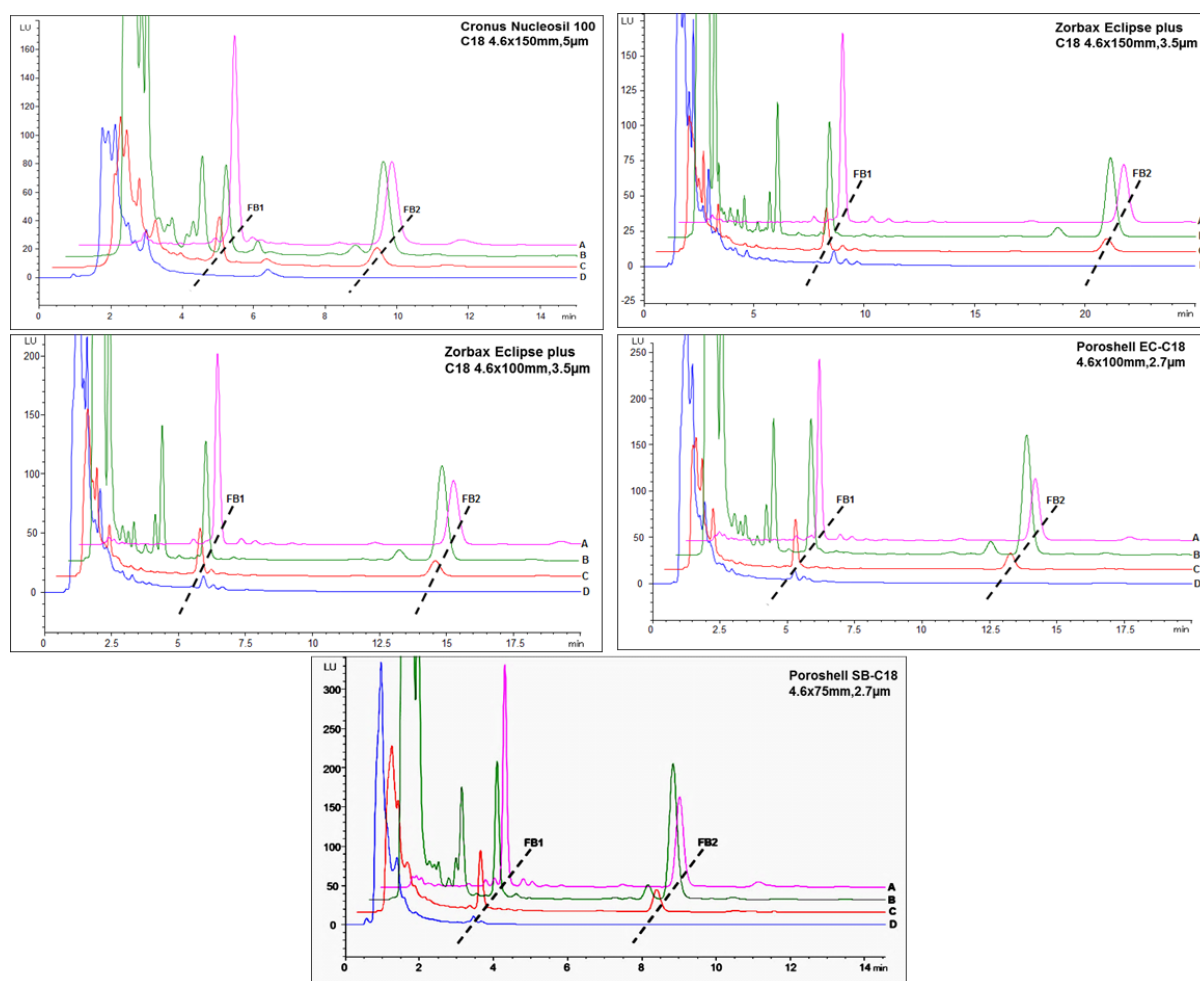

**Figure 3.** HPLC-FLD chromatograms of different matrices contaminated with fumonisins B<sub>1</sub> and B<sub>2</sub> using all tested columns at wavelengths  $\lambda_{ex}$  335nm and  $\lambda_{em}$  440nm (A) mix of fumonisins B<sub>1</sub> and B<sub>2</sub> standard ; (B) fumonisins B<sub>1</sub> and B<sub>2</sub> produced by *Fusarium verticillioides* grown on FIM; (C) fumonisins B<sub>1</sub> and B<sub>2</sub> in spiked cornflakes sample; (D) fumonisins B<sub>1</sub> and B<sub>2</sub> in control cornflakes sample (Non Detected).



The shortest retention times were obtained using the new Agilent Poroshell 120 (75 mm, 2.7 $\mu$ m). Retention times were 3.2 and 7.7 min for both FB<sub>1</sub> and FB<sub>2</sub>, respectively. In contrast, Poroshell columns showed the highest backpressure when compared with Nucleosil and Zorbax eclipse columns.

Generally, the peak width of FB<sub>2</sub> was around twice wider than that of FB<sub>1</sub> using all selected C<sub>18</sub> columns. Poroshell columns, notably the shorter one (75mm) showed the narrowest peaks of both FB<sub>1</sub> and FB<sub>2</sub> (0.12 and 0.24 min), referring to how sharp the peaks were. The reductions in particle size (Cronus Nucleosil 100, 5  $\mu$ m to Zorbax Eclipse plus, 3.5  $\mu$ m) led to 114% and 156% increase in N of FB<sub>1</sub> and FB<sub>2</sub> respectively and subsequently 53% and 61% reduction in the HEPT. Also, the length increases led to an increase in the number of theoretical plates at the same particle sizes (N)(Table 1).

Results for the reduced plate height (h) are also shown in Table 1. The highest h value was for Cronus Nucleosil 100, 5  $\mu$ m recording >3.2 for both FB<sub>1</sub> and FB<sub>2</sub>. However, Zorbax and Poroshell columns had close h values with only differences <0.25. Finally, the shortest time of whole analysis for either cornflakes or culture samples was obtained using both Cronus 150 mm and Agilent 75 mm (105 and 90 min, respectively).

### 3.3. Linearity and sensitivity.

Table 2 illustrates the linear regression equations and the coefficient of determination (R<sup>2</sup>) of FB<sub>1</sub> and FB<sub>2</sub> standards analyzed by different columns. In general, the slope of FB<sub>1</sub> (16.5-17.5) was higher than that of FB<sub>2</sub> (11-13.2) for all columns, referring to the higher FLD response. The coefficient of determination (R<sup>2</sup>) of both fumonisin types using the five columns was very near from 1 (>0.997), indicating an excellent linear response of the FLD detector.

**Table 2.** Linear regression of FB<sub>1</sub> and FB<sub>2</sub> standards using FLD detector and different HPLC C<sub>18</sub> columns.

|                 |                | Cronus<br>150x4.6mm-<br>5 $\mu$ m | Zorbax<br>150x4.6mm-<br>3.5 $\mu$ m | Zorbax<br>100x4.6mm-<br>3.5 $\mu$ m | Poroshell<br>100x4.6mm-<br>2.7 $\mu$ m | Poroshell<br>75x4.6mm-2.7 $\mu$ m |
|-----------------|----------------|-----------------------------------|-------------------------------------|-------------------------------------|--|-----------------------------------|
| FB <sub>1</sub> | Equation       | y = 16.634x - 6.5329              | y = 16.524x - 16.345                | y = 17.489x - 16.421                | y = 16.916x - 22.72                    | y = 17.016x - 5.0586              |
|                 | R <sup>2</sup> | 0.9996                            | 0.9988                              | 0.9981                              | 0.9983                                 | 0.9997                            |
| FB <sub>2</sub> | Equation       | y = 12.46x - 5.3169               | y = 11.079x - 19.252                | y = 12.577x - 6.6159                | y = 12.91x - 22.395                    | y = 13.173x - 19.07               |
|                 | R <sup>2</sup> | 0.9998                            | 0.9966                              | 0.9996                              | 0.9982                                 | 0.9981                            |

### 3.4. Limit of detection and quantification.

LOD and LOQ of FB<sub>1</sub> and FB<sub>2</sub> in both standards solution ( $\mu$ g l<sup>-1</sup>) and sample matrices ( $\mu$ g kg<sup>-1</sup>) were calculated according to Miller and Miller [43], and the data are illustrated in Table 3. Generally, both LOD and LOQ values of FB<sub>2</sub> were more than twice times higher than those of FB<sub>1</sub>, referring to the lower response by FLD detector. The lowest values of both parameters were observed in the case of standards followed by cornflakes and culture. Poroshell 75x4.6mm-2.7 $\mu$ m was the most sensitive column notably with FB<sub>2</sub> type recording LOD<1.2  $\mu$ g l<sup>-1</sup> for standard and <8.5 for culture samples. The same trend was observed regarding to the LOQ parameter. On the other hand, Cronus 150x4.6mm-5 $\mu$ m recorded the highest LOD and LQD of FB<sub>1</sub> and FB<sub>2</sub> in standard. Otherwise, these limits were very close in other substances using the tested columns.

**Table 3.** Comparison of the LOD and LOQ ( $\mu\text{g kg}^{-1}$ ) of Fumonisin B<sub>1</sub> and B<sub>2</sub> (FB<sub>1</sub>, FB<sub>2</sub>) obtained with the columns assayed.

|     |                 |                                      | Cronus<br>150x4.6mm-<br>5 $\mu\text{m}$ | Zorbax<br>150x4.6mm-<br>3.5 $\mu\text{m}$ | Zorbax<br>100x4.6mm-<br>3.5 $\mu\text{m}$ | Poroshell<br>100x4.6mm-<br>2.7 $\mu\text{m}$ | Poroshell<br>75x4.6mm-<br>2.7 $\mu\text{m}$ |
|-----|-----------------|--------------------------------------|---|---|---|--|---|
| LOD | FB <sub>1</sub> | Std. ( $\mu\text{g l}^{-1}$ )        | 1.172                                   | 0.53                                      | 0.94                                      | 0.74   | 0.49  |
|     |                 | Cornflakes ( $\mu\text{g kg}^{-1}$ ) | 1.90                                    | 1.49                                      | 2.64                                      | 4.10   | 1.99  |
|     |                 | Culture ( $\mu\text{g kg}^{-1}$ )    | 3.87                                    | 3.55                                      | 4.66                                      | 4.79   | 3.41  |
|     | FB <sub>2</sub> | Std. ( $\mu\text{g l}^{-1}$ )        | 2.99                                    | 1.73                                      | 2.81                                      | 2.10   | 1.180                                       |
|     |                 | Cornflakes ( $\mu\text{g kg}^{-1}$ ) | 4.83                                    | 4.88                                      | 7.89                                      | 11.70  | 4.85  |
|     |                 | Culture ( $\mu\text{g kg}^{-1}$ )    | 9.87                                    | 11.64                                     | 13.89                                     | 13.65  | 8.43  |
| LOQ | FB <sub>1</sub> | Std. ( $\mu\text{g l}^{-1}$ )        | 3.91                                    | 1.76                                      | 3.15                                      | 2.46   | 1.62  |
|     |                 | Cornflakes ( $\mu\text{g kg}^{-1}$ ) | 6.32                                    | 4.96                                      | 8.8                                       | 13.67  | 6.64  |
|     |                 | Culture ( $\mu\text{g kg}^{-1}$ )    | 12.94                                   | 11.83                                     | 15.52                                     | 15.96  | 11.54                                       |
|     | FB <sub>2</sub> | Std. ( $\mu\text{g l}^{-1}$ )        | 9.97                                    | 5.76                                      | 9.40                                      | 7.00   | 3.93  |
|     |                 | Cornflakes ( $\mu\text{g kg}^{-1}$ ) | 16.10                                   | 16.27                                     | 26.26                                     | 39.00  | 16.16                                       |
|     |                 | Culture ( $\mu\text{g kg}^{-1}$ )    | 32.89                                   | 38.79                                     | 46.28                                     | 45.51  | 28.09                                       |

### 3.5. Analysis of FB<sub>1</sub> and FB<sub>2</sub> in fungal cultures and cornflakes samples.

Table 4 shows fumonisins levels in cultures of 2 isolates of *F. verticilloides* and in cornflakes samples separated by different columns. In general, the isolate S40 was able to produce a high amount of FB<sub>1</sub> and FB<sub>2</sub> ( $>27$  and  $>64 \mu\text{g g}^{-1}$  media, respectively) whereas, a small amount of FB<sub>1</sub> ( $<2.6 \mu\text{g g}^{-1}$  media) and FB<sub>2</sub> ( $<1.3 \mu\text{g g}^{-1}$  media) were produced by the isolate S30. The records of both FB<sub>1</sub> and FB<sub>2</sub> levels in S30 culture had no significant Cronus 150x4.6mm-5 $\mu\text{m}$  was statistically lower than that of Poroshell 100 at  $P < 0.05$ . Control samples of cornflakes contained only low levels of FB<sub>1</sub> ( $<0.14\text{ppm}$ ), which varied among the tested columns.

**Table 4.** Fumonisin B<sub>1</sub> and B<sub>2</sub> levels in fumonisins producing cultures and cornflakes samples.

| C18 Column type                       |                 | $\mu\text{g g}^{-1}$ media $\pm$ SE         |   |                                | Spiked cornflakes (mg kg <sup>-1</sup> $\pm$ SE) |                     |                               |                     |                               |                     |           |
|---------------------------------------|-----------------|---|---|--------------------------------|--|---------------------|-------------------------------|---------------------|-------------------------------|---------------------|-----------|
|                                       |                 | <i>F. verticillioide</i><br><i>s</i><br>S30 | <i>F. verticillioide</i><br><i>s</i><br>S40 |                                | Control  | 0.625               |                               | 1.25                |                               | 2.5                 |           |
|                                       |                 |   |   |                                |  | Calculated<br>Conc. | Recovery%                     | Calculated<br>Conc. | Recovery%                     | Calculated<br>Conc. | Recovery% |
| Cronus 150x4.6mm-5 $\mu\text{m}$      | FB <sub>1</sub> | 2.12 <sup>a</sup> $\pm$ 0.425               | 26.9 <sup>b</sup> $\pm$ 2.50                | ND                             | 0.67 <sup>b</sup> $\pm$ 0.06                     | 107.2               | 1.13 <sup>c</sup> $\pm$ 0.07  | 90.4                | 2.19 <sup>c</sup> $\pm$ 0.10  | 87.6                |           |
|                                       | FB <sub>2</sub> | 0.94 <sup>a</sup> $\pm$ 0.219               | 64.38 <sup>a</sup> $\pm$ 2.81               | ND                             | 0.71 <sup>b</sup> $\pm$ 0.07                     | 113.6               | 1.24 <sup>a</sup> $\pm$ 0.05  | 99.2                | 2.40 <sup>b</sup> $\pm$ 0.11  | 96.2                |           |
| Zorbax 150x4.6mm-3.5 $\mu\text{m}$    | FB <sub>1</sub> | 2.27 <sup>a</sup> $\pm$ 0.45                | 30.51 <sup>ab</sup> $\pm$ 2.04              | 0.032 <sup>b</sup> $\pm$ 0.005 | 0.73 <sup>b</sup> $\pm$ 0.07                     | 112.0               | 1.33 <sup>b</sup> $\pm$ 0.08  | 104.0               | 2.42 <sup>ab</sup> $\pm$ 0.12 | 95.6                |           |
|                                       | FB <sub>2</sub> | 1.17 <sup>a</sup> $\pm$ 0.133               | 67.07 <sup>a</sup> $\pm$ 2.95               | ND                             | 0.68 <sup>a</sup> $\pm$ 0.80                     | 108.8               | 1.26 <sup>a</sup> $\pm$ 0.13  | 100.8               | 2.51 <sup>ab</sup> $\pm$ 0.11 | 100.4               |           |
| Zorbax 100x4.6mm-3.5 $\mu\text{m}$    | FB <sub>1</sub> | 2.57 <sup>a</sup> $\pm$ 0.51                | 30.35 <sup>ab</sup> $\pm$ 1.64              | 0.038 <sup>b</sup> $\pm$ 0.007 | 0.56 <sup>c</sup> $\pm$ 0.07                     | 84.0                | 1.36 <sup>b</sup> $\pm$ 0.10  | 105.8               | 2.34 <sup>b</sup> $\pm$ 0.09  | 93.6                |           |
|                                       | FB <sub>2</sub> | 1.29 <sup>a</sup> $\pm$ 0.32                | 71.58 <sup>a</sup> $\pm$ 3.84               | ND                             | 0.50 <sup>c</sup> $\pm$ 0.03                     | 80.0                | 1.33 <sup>a</sup> $\pm$ 0.07  | 106.4               | 2.70 <sup>a</sup> $\pm$ 0.18  | 108                 |           |
| Poroshell 100x4.6mm-2.7 $\mu\text{m}$ | FB <sub>1</sub> | 2.55 <sup>a</sup> $\pm$ 0.46                | 35.6 <sup>a</sup> $\pm$ 1.85                | 0.021 <sup>c</sup> $\pm$ 0.003 | 0.75 <sup>ab</sup> $\pm$ 0.10                    | 116.6               | 1.49 <sup>ab</sup> $\pm$ 0.08 | 117.5               | 2.56 <sup>a</sup> $\pm$ 0.06  | 101.6               |           |
|                                       | FB <sub>2</sub> | 1.26 <sup>a</sup> $\pm$ 0.28                | 72.6 <sup>a</sup> $\pm$ 0.40                | ND                             | 0.74 <sup>ab</sup> $\pm$ 0.04                    | 118.4               | 1.28 <sup>a</sup> $\pm$ 0.03  | 102.4               | 2.36 <sup>b</sup> $\pm$ 0.05  | 94.4                |           |
| Poroshell 75x4.6mm-2.7 $\mu\text{m}$  | FB <sub>1</sub> | 2.47 <sup>a</sup> $\pm$ 0.51                | 30.74 <sup>ab</sup> $\pm$ 1.67              | 0.140 <sup>a</sup> $\pm$ 0.04  | 0.77 <sup>a</sup> $\pm$ 0.13                     | 101.0               | 1.57 <sup>a</sup> $\pm$ 0.12  | 114.8               | 2.55 <sup>a</sup> $\pm$ 0.15  | 96.5                |           |
|                                       | FB <sub>2</sub> | 1.20 <sup>a</sup> $\pm$ 0.22                | 65.51 <sup>a</sup> $\pm$ 2.12               | ND                             | 0.73 <sup>b</sup> $\pm$ 0.06                     | 116.8               | 1.30 <sup>a</sup> $\pm$ 0.12  | 104                 | 2.56 <sup>ab</sup> $\pm$ 0.09 | 102.4               |           |

Means followed by different superscript letters of the same toxin type within columns are statistically different at  $p < 0.05$ . \*: Recovery= (Actual concentration - Control)/prepared concentration x100.



There were significant variations in fumonisins concentrations of spiked samples. In general, the recoveries of both fumonisins types using all columns were <120%. Both Poroshell columns recorded the highest recoveries, notably at 0.625 and 1.25 levels, whereas close values to 100% were obtained at 2.5 mg kg<sup>-1</sup>.

### 3.6. Discussion

There are no available or similar studies comparing the suitability of different chromatographic analytical columns for both fumonisins and have compared the analysis performance using different matrices. In our aim to provide optimal analysis conditions for many laboratories around the world where access to UPLC or LC/MS-MS technologies are just unaffordable and in continuation with our previous work on ochratoxin A [40], the current study has considered the suitability of 5 different C<sub>18</sub> columns with different particle sizes and porosities to analyze FB<sub>1</sub> and FB<sub>2</sub>. In addition to standard solutions, fungal cultures of fumonisins-producing fungi and cornflakes were chosen for column evaluation to provide evidence of how the selected methodology can be applied in a wide range of applications spanning from mycological research to the analysis of food matrices. An extensive chromatographic comparison has been carried out, exploring the separation of both FB<sub>1</sub> and FB<sub>2</sub>.

The derivatization step, which needs to be performed just before each injection, was one of the challenges we faced in this study, as it hampered the potential use of autosampler devices. However, after our initial literature review, it became clear that the use of OPA was the best option in terms of LOD and LOQ [33]. Finally, OPA was selected in an attempt to save analysis time and solvents without affecting the quality of separation. Thus, this is a step where future research needs to be further developed and could be a good improvement target for analytical laboratories and research institutions.

Although the retention times of both FB<sub>1</sub> and FB<sub>2</sub> using the modified mobile phase were higher than those using AOAC mobile phase, a better resolution was observed (Figure 2). Increasing the aqueous proportion from 23% to 25% at the expense of methanol proportion increased the polarity of the mobile phase. Consequently, fumonisins were held on the column for a longer time.

The records of the tested columns' backpressure were in line with their properties. It decreased in the shorter column of the same particle size either within Zorbax or Poroshell, whereas smaller particle sizes of the same column length also had higher backpressure (150mm: from 5 to 3.5 µm).

Carbon load and surface area played an important role in the retention times of FB<sub>1</sub> and FB<sub>2</sub>. Albeit Cronus 150 mm had the biggest particle sizes (5µm) among the examined columns, the separation was faster than that obtained by both Zorbax Columns (150 and 100mm:3.5µm). Silica particles of the Cronus column are loaded with 14% carbon and occupy 350m<sup>2</sup> g<sup>-1</sup> surface area exceeding those of Zorbax Eclipse Plus columns (9% and 160 m<sup>2</sup> g<sup>-1</sup>). This made FB<sub>1</sub> and FB<sub>2</sub> passed through Cronus particles quicker than through the Zorbax bed. Separation using Cronus was even faster than the new Poroshell 120 EC-C<sub>18</sub> 4.6x100mm, 2.7µm column packed with core fused particles (0.5µm thick of porous outer layer and 1.7µm solid-core silica). This is because of the lower carbon load and surface area of the latter column, 8%, and 120 m<sup>2</sup>g<sup>-1</sup>, respectively [45]. Only the shorter poroshell column, 75 mm, was faster than Cronus in fumonisins separation (3.2 and 7.7 min for FB<sub>1</sub> and FB<sub>2</sub>, respectively).

Similar to Cronus, Discovery C18 150x4.5 mm, 5  $\mu\text{m}$  column was used by Solfrizzo *et al.* [37] to separate FB<sub>1</sub> and FB<sub>2</sub> after cleaning up the cornflakes samples IAC cartridge at 7 and 17.1 min, respectively. They were also eluted after 6.5 and 11.5 min using ZORBAX Eclipse® XDB C18 column 150x4.5mm,  $\mu\text{m}$  [12].

Ndube *et al.* [33] separated FB<sub>1</sub> and FB<sub>2</sub> from maize samples after derivatization with naphthalene-2,3-dicarboxaldehyde and SAX column clean up. Although they used Phenomenex, Luna C18 5  $\mu\text{m}$  column (75 mm  $\times$  4.60 mm), which was half-length of Cronus column, more retention time was required to elute FB<sub>1</sub> and FB<sub>2</sub> (7.2 and 17.5min).

The peak width of both FB<sub>1</sub> and FB<sub>2</sub> reflected how the sharpness of these peaks was. Poroshell columns gave the narrowest width, notably the shorter one (75mm). This is because the fumonisins molecules move a short distance through the fused core particles when compared with the completely porous particles.

Both retention time and peak width can be used to identify the number of theoretical plates, which was the highest using Zorbax Eclipse 150mm 3.5 mm. It had the broadest width and the longest retention time. Consequently, N numbers were used to calculate HEPT based on the column length, which was the highest in the Cronus column. However, none of these parameters can be used to evaluate the efficiency of the tested column alone or together. The only reduced plate height (h) parameter can be used to evaluate as it considers the column particle size. Cronus had the highest h value referring to the well-packing process with the ideal packed bed. This can be interpreted by the fact that the smaller the particle size is, the greater the difficulty in preparing a well-packed column bed is [46]. This rule is not applicable in the case of different porosity columns with different lengths like Zorbax and Poroshell columns, which had close h values for FB<sub>1</sub> and FB<sub>2</sub>.

Linearity expressed by  $R^2$  in determination for FB<sub>1</sub> and FB<sub>2</sub> is the first evaluation step of the column performance. All tested columns exhibited good linearity in addition to the small intercept values. A similar finding was observed, recording higher than 0.9997 for both toxins [12].

LOD and LOQ are important parameters for any eluents or mycotoxins like fumonisins that contaminate corn and corn products at high levels [47]. LOQ should be less than the permissible limits of FB<sub>1</sub> and FB<sub>2</sub> that reach 200  $\mu\text{g kg}^{-1}$  for babies' food [26]. The extraction method and matrix type noticeably affected LOQ using all tested columns. The lowest values were for the standard solution followed by cornflakes samples which were cleaned up by C18 cartridge, then by culture samples. LOQ for both toxins in cornflakes samples was  $< 23 \mu\text{g kg}^{-1}$  using Cronus, Zorbax 150, and Poroshell 75 columns, whereas it recorded 35 and 53  $\mu\text{g kg}^{-1}$  using Zorbax 100 and Poroshell 100, respectively. A close value was obtained by De Girolamo *et al.* [48], who estimated LOQ for both FB<sub>1</sub> and FB<sub>2</sub> in masa flour by 25  $\mu\text{g kg}^{-1}$  using Symmetry Shield (150 mm $\times$ 4.6 mm, 5  $\mu\text{m}$ ) column. In contrast, it was 400  $\mu\text{g kg}^{-1}$  in milled corn using IAC for cleaning up and MS detector [1]. In general, all values in the present study are much lower than the permissible limit of fumonisins in cornflakes, 800  $\mu\text{g kg}^{-1}$  [26]. Concerning the cornflake sample, it was naturally contaminated with FB<sub>1</sub>, recording  $< 140 \mu\text{g kg}^{-1}$  which was much lower than the permissible limits established by EC [26]. That value is very close to that obtained by Solfrizzo *et al.* [37], who calculated the contamination levels of both FB<sub>1</sub> and FB<sub>2</sub> in 18 cornflake samples in the average of 0.157 and 0.036  $\mu\text{g kg}^{-1}$ , respectively. Although all columns achieved good performance regarding the LOQ ( $< 140 \mu\text{g kg}^{-1}$  for FB<sub>1</sub>), there was a variation in the detection of the FB<sub>1</sub> between the tested column. A

clear peak of FB<sub>1</sub> was observed using Poroshell 75 column, whereas others had an unclear peak, notably using the Cronus column (Figure 3).

In general, there was a small up-shift in the obtained recoveries (<8%) than the acceptance criteria, 70%-110%, established by the European Commission [49]. This can be attributed to the matrix effect of the control cornflake samples, which chelated toxins. However, the spiking process can stimulate releasing the chelated toxin. Contrarily, Solfrizzo *et al.* [37] mentioned that low recoveries of FB<sub>1</sub> and FB<sub>2</sub> in corn-based foods as cornflakes (<30%) were obtained using the SAX column as a clean-up procedure. They explained that the existing iron in commercial cornflakes could chelate the free fumonisin at the two tricarballic acid groups, which would not be able to hold in the SAX ion exchange column.

Except for Zorbax 100, spiked cornflake samples' recoveries at 0.625 mg kg<sup>-1</sup> were higher than those at 1.25 and 2.5 mg kg<sup>-1</sup>. This finding agrees with Li *et al.* [14], who reported that the recovery levels of spiked corn samples decreased by increasing the spiking levels. A similar study achieved 102.6% and 95.1% recoveries for both FB<sub>1</sub> and FB<sub>2</sub> in spiked cornflakes samples [37]. Likewise, Muscarella *et al.* [12] found that the recovery values ranged from 87 to 94% for FB<sub>1</sub> and 70 to 75% for FB<sub>2</sub> in cornflake samples.

#### 4. Conclusions

Using C<sub>18</sub> cartridge in cleaning up followed by OPA derivatization was an effective method in FB<sub>1</sub> and FB<sub>2</sub> determination in cornflakes. Albeit Nucleosil Cronus 150mmx4.6mm, 5 µm had porous particles. It had the closest results to Poroshell 75 in saving time, solvents and gave good recovery value with relatively low pressure. Meanwhile, separating FB<sub>1</sub> and FB<sub>2</sub> using solid core particle techniques in poroshell columns did not make a big difference compared with end-capping Nucleosil Cronus. Finally, column diameter and its particle type and size play an important role in fumonisin B<sub>1</sub> and B<sub>2</sub> analysis.

#### Funding

This research received no external funding.

#### Acknowledgments

The authors thank Dr. Sejakosi Mohale for providing the fungal strains.

#### Conflicts of Interest

The authors declare no conflict of interest.

#### References

1. Dohnal, V.; Ježková, A.; Polišenská, I.; Kuca, K. Determination of Fumonisin in Milled Corn Grains Using HPLC–MS. *Journal of Chromatographic Science* **2010**, *48*, 680–684, <https://doi.org/10.1093/chromsci/48.8.680>.
2. Abd-El Fatah, S.I.; Naguib, M.M.; El-Hossiny, E.N.; Sultan, Y.Y.; Abodalam, T.H.; Yli-Mattila, T. Molecular versus Morphological Identification of *Fusarium spp.* isolated from Egyptian corn. *Research Journal of Pharmaceutical, Biological and Chemical Sciences* **2015**, *6*, 1813–1822, [http://rjpbcs.com/.../\[243\].pdf](http://rjpbcs.com/.../[243].pdf).
3. Beccari, G.; Stępień, Ł.; Onofri, A.; Lattanzio, V.M.T.; Ciasca, B.; Abd-El Fatah, S.I.; Valente, F.; Urbaniak, M.; Covarelli, L. *In Vitro* fumonisin biosynthesis and genetic structure of *Fusarium verticillioides* strains

- from five Mediterranean countries. *Microorganisms* **2020**, *8*, 241, <https://doi.org/10.3390/microorganisms8020241>.
4. Ibrahim, N.F.; Mohd, M.H.; Nor, N.M.I.M.; Zakaria, L. Mycotoxigenic potential of *Fusarium* species associated with pineapple diseases. *Archives of Phytopathology and Plant Protection* **2020**, *53*, 216–228, <https://doi.org/10.1080/03235408.2020.1736971>.
5. Chavez, R.A.; Cheng, X.; Stasiewicz, M.J. A review of the methodology of analyzing aflatoxin and fumonisin in single corn kernels and the potential impacts of these methods on food security. *Foods* **2020**, *9*, 297, <https://doi.org/10.3390/foods9030297>.
6. IARC. Fumonisin B1. In: International Agency for Research on Cancer monographs on the evaluation of the carcinogenic risks to humans: Some traditional herbal medicines, some mycotoxins, naphthalene and styrene. IARC, Lyon, France, **2002**, *82*, 301–366.
7. Yu S., Jia B., Liu N., Yu D., Wu A. Evaluation of the individual and combined toxicity of fumonisin mycotoxins in human gastric epithelial cells. *International Journal of Molecular Sciences* **2020**, *21*, 5917–5930, <https://doi.org/10.3390/ijms21165917>.
8. Voss, K.; Riley, R. Fumonisin Toxicity and Mechanism of Action: Overview and Current Perspectives. *Food Safety* **2013**, *1*, 49–69, <https://doi.org/10.14252/foodsafetyfscj.2013006>.
9. Rao, Z.-X.; Tokach, M.D.; Woodworth, J.C.; DeRouchey, J.M.; Goodband, R.D.; Calderón, H.I.; Dritz, S.S. Effects of Fumonisin-Contaminated Corn on Growth Performance of 9 to 28 kg Nursery Pigs. *Toxins* **2020**, *12*, 604, <https://doi.org/10.3390/toxins12090604>.
10. Arumugam, T.; Ghazi, T.; Chuturgoon, A.A. Molecular and epigenetic modes of Fumonisin B<sub>1</sub> mediated toxicity and carcinogenesis and detoxification strategies. *Critical Reviews in Toxicology* **2021**, *51*, 76–94, <https://doi.org/10.1080/10408444.2021.1881040>.
11. Yayeh, T.; Jeong, H.R.; Park, Y.S.; Moon, S.; Sur, B.; Yoo, H.-S.; Oh, S. Fumonisin B1-Induced Toxicity Was Not Exacerbated in Glutathione Peroxidase-1/Catalase Double Knock Out Mice. *Biomolecules and Therapeutics* **2021**, *29*, 52–57, <https://doi.org/10.4062/biomolther.2020.062>.
12. Muscarella, M.; Lo Magro, S.; Nardiello, D.; Palermo, C.; Centonze, D. Development of a new analytical method for the determination of fumonisins B1 and B2 in food products based on high performance liquid chromatography and fluorimetric detection with post-column derivatization. *Journal of Chromatography A* **2008**, *1203*, 88–93, <https://doi.org/10.1016/j.chroma.2008.07.034>.
13. Chen, J.; Wei, Z.; Wang, Y.; Long, M.; Wu, W.; Kuca, K. Fumonisin B1: Mechanisms of toxicity and biological detoxification progress in animals. *Food and Chemical Toxicology* **2021**, *149*, 111977, <https://doi.org/10.1016/j.fct.2021.111977>.
14. Li, R.; Tao, B.; Pang, M.; Liu, Y.; Dong, J. Natural occurrence of fumonisins B1 and B2 in maize from three main maize-producing provinces in China. *Food Control* **2015**, *50*, 838–842, <https://doi.org/10.1016/j.foodcont.2014.09.034>.
15. García-Díaz, M.; Gil-Serna, J.; Vázquez, C.; Botia, M.N.; Patiño, B. A comprehensive study on the occurrence of mycotoxins and their producing fungi during the Maize production cycle in Spain. *Microorganisms* **2020**, *8*, 141, <https://doi.org/10.3390/microorganisms8010141>.
16. Niyibituronsa, M.; Mukantwali, C.; Nzamwita, M.; Hagenimana, G.; Niyoyita, S.; Niyonshima, A.; Hakizimana, C.; Ndilu, L.; Nyirahanganyamunsi, G.; Nkurunziza, E.; Sendegeya, P.; Niyonteze, G.; Muhutu, J.C.; Shingiro, J.B.; Umuhire, J.; Nyirahorana, C.; Ingabire, A.C.; Nyiranshuti, A.; Sibomana, E.; Usabyembabazi, M.; Munganyinka, E.; Gasana, A.P.; Kamaraba, I.; Izamuhaye, J.C.; Nkundanyirazo, E.; Ngaboyisonga, C.; Bucagu, C.; Karangwa, P. Assessment of aflatoxin and fumonisin contamination levels in maize and mycotoxins awareness and risk factors in Rwanda. *African Journal of Food, Agriculture, Nutrition and Development* **2020**, *20*, 16420–16446, <https://doi.org/10.18697/ajfand.93.19460>.
17. Alves dos Reis, T.; Zorzete, P.; Pozzi, C.R.; Nascimento da Silva, V.; Ortega, E.; Corrêa, B. Mycoflora and fumonisin contamination in Brazilian sorghum from sowing to harvest. *Journal of the Science of Food and Agriculture* **2010**, *90*, 1445–1451, <https://doi.org/10.1002/jsfa.3962>.
18. Trucksess, M.W.; Cho, T.H.; Ready, D.E. Liquid chromatographic method for fumonisin B1 in sorghum syrup and corn-based breakfast cereals. *Food additives and Contaminants* **2000**, *17*, 161–166, <https://doi.org/10.1080/026520300283513>.
19. Tseng, T.C.; Tu, J.C. Mycoflora and mycotoxins in adzuki and mung beans produced in Ontario, Canada. *Micobios*, **1997**, *90*, 87–95.
20. Mallmann, C.A.; Santurio, J.M.; Almeida, C.A.A.; Dilkin, P. Fumonisin B1 levels in cereals and feeds from southern Brazil. *Arquivos do Instituto Biológico* **2001**, *68*, 41–45.

21. Iqbal, S.Z.; Rehman, B.; Jinap, S.; Akram, N.; Ahmad, M.N.; Sanny, M. Assessment of fumonisin B<sub>1</sub> concentrations in wheat and barley products in the Punjab region of Pakistan. *Journal of Food Protection* **2020**, *83*, 1284-1288, <https://doi.org/10.4315/0362-028X.JFP-19-361>.
22. Martins, M.L.; Martins, H.M.; Bernardo, F. Fumonisin B<sub>1</sub> and B<sub>2</sub> in black tea and medicinal plants. *Journal of Food Protection* **2001**, *64*, 1268-1270, <https://doi.org/10.4315/0362-028X-64.8.1268>.
23. Ikeagwulonu, R. C., C. Onyenekwe, C., O. Oshim, I., Ikeagwulonu, Z. C., K. Ojide, C., & A. Olise, N. Evaluation of the Levels of Ochratoxin A and Fumonisin B<sub>1</sub> Mycotoxins in Herbal Traditional Medicines Selected from Vendor Dealers in Ebonyi State of Nigeria. *Journal of Advances in Medical and Pharmaceutical Sciences* **2020**, *22*, 1-5, <https://doi.org/10.9734/jamps/2020/v22i130149>.
24. Becker-Algeri, T.A.; Heidtmann-Bemvenuti, R.; Hackbart, H.C.; Badiale-Furlong, E. Thermal treatments and their effects on the fumonisin B level in rice. *Food Control* **2010**, *2*, 488-493, <https://doi.org/10.1016/j.foodcont.2013.05.016>.
25. CAST. Mycotoxins: risks in plant animal and human systems. Ames, Iowa: Council for Agricultural Science and Technology, **2003**, 199.
26. EC. European Commission, Commission Directive 2007, 1126/2007 EC, Amending Regulation (EC) No 1881/2006, **2007**.
27. Belajová, E.; Rauová, D. Single laboratory-validated HPLC methods for determination of ochratoxin A, fumonisin B<sub>1</sub> and B<sub>2</sub>, zearalenone and deoxynivalenol in cereals and cereal-based foods, *Journal of Food and Nutrition Research* **2010**, *49*, 57-68.
28. Silva, L.; Fernández-Franzón, M.; Font, G.; Pena, A.; Silveira, I.; Lino, C.; Mañes, J. Analysis of fumonisins in corn-based food by liquid chromatography with fluorescence and mass spectrometry detectors. *Food Chemistry* **2009**, *112*, 1031-1037, <https://doi.org/10.1016/j.foodchem.2008.06.080>.
29. Jo, H.W.; Park, M.K.; Heo, H.m.; Jeon, H.J.; Choi, S.D.; Lee, S.E.; Moon J.K. Simultaneous determination of 13 mycotoxins in feedstuffs using QuEChERS extraction. *Applied Biological Chemistry* **2021**, *64*, <https://doi.org/10.1186/s13765-021-00602-9>.
30. Nakagawa, H.; Hashimoto, R.; Matsuo, Y.; Sago, Y.; Yokoyama, K.; Takahashi, H. Detection and Determination of Fumonisin B<sub>1</sub>, B<sub>2</sub>, and B<sub>3</sub> Contaminating Japanese Domestic Wine by Liquid Chromatography Coupled to Tandem Mass Spectrometry (LC-MS/MS). *Current Microbiology* **2020**, *77*, 3057-3064, <https://doi.org/10.1007/s00284-020-02113-0>.
31. Smaoui, S.; Ben Braïek O.; Ben Hlima, H. Mycotoxins Analysis in Cereals and Related Foodstuffs by Liquid Chromatography-Tandem Mass Spectrometry Techniques. *Journal of Food Quality* **2020**, <https://doi.org/10.1155/2020/8888117>.
32. Greeff-Laubscher, M.R.; Beukes, I.; Marais, G.J.; Jacobs, K. Mycotoxin production by three different toxigenic fungi genera on formulated abalone feed and the effect of an aquatic environment on fumonisins. *Mycology* **2020**, *11*, 105-117, <https://doi.org/10.1080/21501203.2019.1604575>.
33. Ndubea, N.; van der Westhuizen, L.; Green, I.R.; Shephard, G.S. HPLC determination of fumonisin mycotoxins in maize: A comparative study of naphthalene-2,3-dicarboxaldehyde and o-phthalaldehyde derivatization reagents for fluorescence and diode array detection. *Journal of Chromatography B* **2001**, *879*, 2239-2243, <https://doi.org/10.1016/j.jchromb.2011.06.007>.
34. Dombrink-Kurtzman, M.A.; Dvorak, T.J. Fumonisin Content in Masa and Tortillas from Mexico. *Journal of Agriculture and Food Chemistry* **1999**, *47*, 622-627, <https://doi.org/10.1021/jf9807162>.
35. Kim, E.K.; Scott, P.M.; Lau, B.P.; Lewis, D.A. Extraction of fumonisins B<sub>1</sub> and B<sub>2</sub> from white rice flour and their stability in white rice flour, cornstarch, cornmeal, and glucose. *Journal of Agriculture and Food Chemistry* **2002**, *50*, 3614-3620, <https://doi.org/10.1021/jf020058s>.
36. AOAC. Natural toxins. In: Official methods of analysis, 21st ed. Association of Official Analytical Chemists, Washington D.C., **2019**.
37. Solfrizzo, M.; De Girolamo, A.; Visconti, A. Determination of fumonisins B<sub>1</sub> and B<sub>2</sub> in cornflakes by high performance liquid chromatography and immunoaffinity clean-up. *Food Additives and Contaminants* **2001**, *18*, 227-235, <https://doi.org/10.1080/02652030010021963>.
38. Riley, R.T.; Torres, O.A.; Palencia, E. International shipping of fumonisins from maize extracts on C18 sorbent. *Food Additives and Contaminants* **2006**, *23*, 826-832, <https://doi.org/10.1080/02652030600699650>.
39. Medina, A.; Magan, N. Comparison of three different C18 HPLC columns with different particle sizes for the optimization of aflatoxins analysis. *Journal of Chromatography B* **2012**, *889-890*, 138-143 <https://doi.org/10.1016/j.jchromb.2012.01.036>.



40. Sultan, Y.; Magan N.; Medina Á. Comparison of five different C18 HPLC analytical columns for the analysis of ochratoxin A in different matrices. *Journal of Chromatography B* **2014**, 971, 89-93, <https://doi.org/10.1016/j.jchromb.2014.09.016>.
41. Jurado, M.; Marín, P.; Magan, N.; González-Jaén, M.T. Relationship between Solute and Matric Potential Stress, Temperature, Growth, and FUM1 Gene Expression in Two *Fusarium verticillioides* Strains from Spain. *Applied and Environmental Microbiology* **2008**, 74, 2032-2036, <https://doi.org/10.1128/AEM.02337-07>.
42. Lazzaro, I.; Falavigna, C.; Dall'Asta, C.; Proctor, R.H.; Galaverna, G.; Battilani, P. Fumonisin B<sub>1</sub>, A and C profile and masking in *Fusarium verticillioides* strains on fumonisin-inducing and maize-based media. *International Journal of Food Microbiology* **2012**, 159, 93-100, <https://doi.org/10.1016/j.ijfoodmicro.2012.08.013>.
43. Miller, J.N.; Miller, J.C. Statistics and Chemometrics for Analytical Chemistry, 4th ed., Prentice Hall, Harlow, **2000**.
44. van Deemter, J.J. Longitudinal diffusion and resistance to mass transfer as causes of nonideality in chromatography. *Chemical Engineering Science* **1995**, 50, 3867, [https://doi.org/10.1016/0009-2509\(96\)81812-4](https://doi.org/10.1016/0009-2509(96)81812-4).
45. Agilent. The LC Handbook Guide to LC Columns and Method Development. Agilent Technologies, Inc. Published in USA. February, 2016, <https://www.agilent.com/cs/library/primers/public/LC-Handbook-Complete-2.pdf>.
46. Feketea, S.; Ganzlera, K.; Fekete, J. Facts and myths about columns packed with sub-3µm and sub-2µm particles. *Journal of Pharmaceutical and Biomedical Analysis* **2010**, 51, 56-64, <https://doi.org/10.1016/j.jpba.2009.08.003>.
47. Shephard, G.S.; Burger, H.; Rheeder, J.P.; Albertsa, J.F.; Gelderblom, W.C.A. The effectiveness of regulatory maximum levels for fumonisin mycotoxins in commercial and subsistence maize crops in South Africa. *Food Control* **2019**, 97, 77-80, <https://doi.org/10.1016/j.foodcont.2018.10.004>.
48. De Girolamo, A.; Pascale, M.; Visconti, A. Comparison of methods and optimisation of the analysis of fumonisins B<sub>1</sub> and B<sub>2</sub> in masa flour, an alkaline cooked corn product. *Food Additives and Contaminants* **2011**, 28, 667-675, <https://doi.org/10.1080/19440049.2011.555846>.
49. EC. European Commission, Commission Regulation (EC) No 401/2006 of 23 February 2006 laying down the methods of sampling and analysis for the official control of the levels of mycotoxins in foodstuffs. *Official Journal of the European Union*, L70, 12e34, **2006**.